

Product Data Sheet

Product Name: PROTAC USP39 Degradar-1

Cat. No.: GC80895

Chemical Properties

Cas. No.


Formula $C_{41}H_{52}ClN_7O_8S_2$ M.Wt 870.48

Solubility DMSO: 100 mg/mL (114.88 mM; Storage Store at -20°C
Need ultrasonic)

General tips For obtaining a higher solubility , please warm the tube at 37 °C and shake it in the ultrasonic bath for a while. Stock solution can be stored below -20°C for several months.

Shipping Condition Evaluation sample solution : ship with blue ice All other available size: ship with RT , or blue ice upon request.

Structure

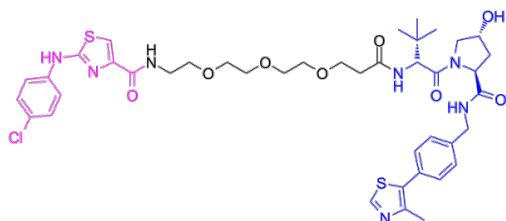

Background

PROTAC USP39 Degradar-1 is a USP39 PROTAC degrader with K_d values of 136 nM and 232 nM,

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

and a DC₅₀ value of 1.06 nM (24 h). PROTAC USP39 Degradator-1 binds to the VHL E3 ubiquitin ligase, forms a ternary complex with USP39, mediates USP39 ubiquitination and proteasomal degradation dependent on VHL recruitment, modulates splicing activity and induces 5' splice site-specific splicing patterns. PROTAC USP39 Degradator-1 can be used in the research of cervical cancer [1]. (Pink: USP39 ligand; Blue: VHL ligand; Black: linker).

In Vitro, PROTAC USP39 Degradator-1 (compound USP39_PROTAC_V1) reduces the thermal stability of purified full-length USP39 protein, decreasing its melting temperature by 1.5 °C[1]. PROTAC USP39 Degradator-1 (0.5 nM-5 µM; 24 h) induces dose-dependent degradation of USP39 protein in HEK293T, LNCap and Kelly cells[1]. PROTAC USP39 Degradator-1 (10 nM-1 µM; 24 h) alters splicing activity in HEK293T cells[1]. PROTAC USP39 Degradator-1 (1 nM-1 µM; 6-24 h) exerts significant target protein degradation activity in HeLa cells transiently overexpressing USP39-OFP, with the optimal treatment condition being 16 h of serum starvation followed by 24 h of compound treatment[1]. PROTAC USP39 Degradator-1 (0.5 nM-5 µM; 6 h-24 h) exerts potent activity in degrading endogenous USP39 protein in HeLa-WT cells (with a 24 h DC₅₀ of approximately 1.06 nM). It degrades its target with high specificity across the whole proteome, resulting in significant downregulation of USP39 and stress-related proteins without affecting homologous proteins[1]. The USP39-degrading effect of PROTAC USP39 Degradator-1 (0.1 nM-0.1 µM; 24 h) is significantly inhibited by the proteasome inhibitor MG132 in co-treated HeLa-WT cells[1]. The degradation effect of PROTAC USP39 Degradator-1 (10 nM; 24 h) is blocked by the Neddylation inhibitor MLN4924 in HeLa-USP39-OFP cells co-treated with 1 µM MLN4924[1]. When PROTAC USP39 Degradator-1 is co-treated with 5 µM MLN4924 in HeLa-USP39-OFP cells, its DC₅₀ increases from 13 nM to 204 nM, and its degradation activity is significantly attenuated[1]. PROTAC USP39 Degradator-1 (0.5 nM-0.5 µM) shows significantly attenuated USP39-degrading activity in VHL-knockout HeLa cells, and exhibits potent degrading activity only in the presence of wild-type VHL in HEK293T cells co-expressing USP39-OFP and different VHL plasmids (VHL_WT or truncated VHL_m)[1]. PROTAC USP39 Degradator-1 (0.5 nM-0.5 µM) promotes the formation of a ternary complex between the target protein and VHL_WT in HEK293T cells transfected with different types of plasmids (KD = 12 nM), whereas the complex formation is extremely weak in the presence of VHL_m[1]. PROTAC USP39 Degradator-1 (100 nM; 6 h) significantly enriches USP39 and polyubiquitination, and promotes the binding of USP39 to VHL in HeLa cells co-treated with 10 µM MLN4924[1]. PROTAC USP39

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

Degrader-1 (10 nM; 24 h) induces aberrant usage of non-canonical 5'-splice sites in HEK293T cells, recapitulating the splicing defects associated with USP39 depletion^[1].

References:

[1]. Schäfer D, et al. Targeting the Spliceosomal Protein USP39 Through Allosteric Ligands and PROTAC-Induced Degradation. *Angewandte Chemie (International ed. in English)*. 2026 Jan 28;65(5):e16809.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA