

Product Data Sheet

Product Name: TPTZ
Cat. No.: GC45581

Chemical Properties

Cas. No. 3682-35-7

SMILES C1(C2=NC(C3=NC=CC=C3)=NC(C4=CC=CC=N4)=N2)=CC=CC=N1


Formula $C_{18}H_{12}N_6$ M.Wt 312.3

Solubility DMF: 10mg/mL, DMSO: 10mg/mL, Ethanol: 20mg/mL, Ethanol:PBS (pH 7.2) (1:4): 0.2mg/mL
Storage Store at -20°C

General tips For obtaining a higher solubility, please warm the tube at 37 °C and shake it in the ultrasonic bath for a while. Stock solution can be stored below -20°C for several months.

Shipping Condition Evaluation sample solution: ship with blue ice. All other available size: ship with RT, or blue ice upon request.

Structure

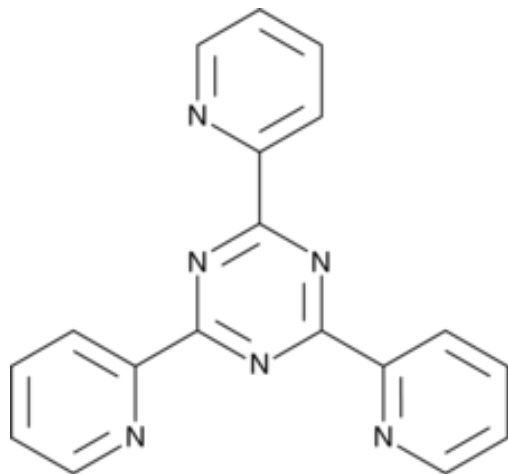

Background

TPTZ is a colorimetric reagent commonly used in the detection of iron. 1,2 TPTZ forms a complex with Fe(II) that can be quantified by colorimetric detection at 594 nm as a measure of iron concentration.

References

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

1. Collins, P.F., and Diehl, H. 2,4,6-Tripyridyl-s-triazine as a reagent for iron. Determination of iron in limestone, silicates, and refractories. Anal. Chem. 31(11), 1862-1867 (1959).
2. Dehghani, Z., Dadfarnia, S., Shabani, A.M.H., et al. Determination of iron species by combination of solvent assisted-dispersive solid phase extraction and spectrophotometry. Anal. Bioanal. Chem. Res. 2(1), 13-21 (2015).

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA