

Product Data Sheet

Product Name: 4-Methylumbelliferyl-β-D-glucuronide hydrate (MUG)

Cat. No.: GC30519

Chemical Properties

Cas. No. 881005-91-0

SMILES O[C@H]([C@H]([C@@H]([C@@H](C(O)=O)O1)O)O)[C@@H]1OC2=CC=C(C(O3)=C2)C(C)=CC3=O.O


Formula $C_{16}H_{18}O_{10}$ M.Wt 370.31

Solubility DMSO : 150 mg/mL (405.07 mM); Water : 5.2 mg/mL (14.04 mM) Storage Store at 4°C, protect from light

General tips For obtaining a higher solubility, please warm the tube at 37 °C and shake it in the ultrasonic bath for a while. Stock solution can be stored below -20°C for several months.

Shipping Condition Evaluation sample solution : ship with blue ice All other available size: ship with RT, or blue ice upon request.

Structure

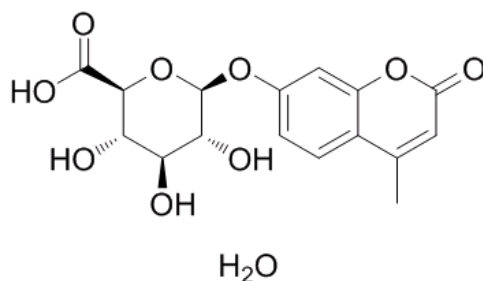

Protocol

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

Kinase experiment: GLUase activity of E. coli cells is measured. Three millilitres of 4-Methylumbelliferyl- β -D-glucuronide hydrate (MUGlu) solution (55 mg of hydrated MUGlu and 20 mL of Triton X-100 in 50 mL of sterile water) is added to each flask (final concentration: 165 mg/L). The incubation temperature is 44°C. One hundred microlitres of 2 M NaOH solution is added to each 2.9-mL aliquot to obtain a pH > 10 before the fluorescence measurement [1].

References:

[1]. George I, et al. Use of beta-D-galactosidase and beta-D-glucuronidase activities for quantitative detection of total and fecal coliforms in wastewater. Can J Microbiol. 2001 Jul;47(7):670-5.

Background

4-Methylumbelliferyl- β -D-glucuronide hydrate is a fluorogenic substrate (λ_{ex} =362 nm, λ_{em} =445 nm).

4-Methylumbelliferyl- β -D-glucuronide hydrate is a fluorogenic substrate with λ_{ex} =362 nm and λ_{em} =445 nm. The half-saturation constants, K_m , that are calculated for sewage samples are close to the K_m values calculated based on experimental plots presented previously for freshwaters [1].

[1]. George I, et al. Use of beta-D-galactosidase and beta-D-glucuronidase activities for quantitative detection of total and fecal coliforms in wastewater. Can J Microbiol. 2001 Jul;47(7):670-5.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA