

Product Data Sheet

Product Name: Hoechst 33258

Cat. No.: GC64380

Chemical Properties

Cas. No. 23491-44-3

Formula $C_{25}H_{24}N_6O$

M.Wt 424.5


Solubility DMSO : 41.67 mg/mL (98.16 mM; ultrasonic and warming and heat to 60°C)

Storage -20°C, protect from light

General tips For obtaining a higher solubility , please warm the tube at 37 °C and shake it in the ultrasonic bath for a while. Stock solution can be stored below -20°C for several months.

Shipping Condition Evaluation sample solution : ship with blue ice All other available size: ship with RT , or blue ice upon request.

Structure

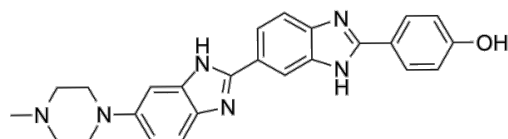

Protocol

This plan only provides a guide, please modify it to meet your specific needs.

1. Prepare Hoechst staining solution

(1) Prepare Hoechst dye stock solution: Use DMSO to dissolve solid and prepare a 10 mg/mL Hoechst dye stock solution.

Note: Hoechst stock solution is recommended to be aliquoted and stored in the dark at -4°C or -20°C to avoid repeated freezing and thawing.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

(2) Working solution preparation: Use preheated serum-free medium or buffer (such as HBSS or PBS) to dilute the stock solution and prepare a Hoechst working solution with a concentration of 10 µg/mL.

Note: Please adjust the concentration of Hoechst working fluid according to the actual situation and prepare it immediately.

2. Cell staining

2.1 Suspension cells (taking 6-well plate as an example)

(1) Centrifuge suspended cells at 1000g for 3-5 minutes. Discard the supernatant and wash twice with PBS for 5 minutes each time.

(2) Add 1 mL of Hoechst dye working solution and incubate at room temperature in the dark for 5-10 minutes.

(3) After the incubation, centrifuge at 1000g for 5 minutes, remove the supernatant, add PBS and wash 2-3 times, 5 minutes each time.

(4) Use serum-free cell culture medium or PBS to resuspend the cells and observe them through fluorescence microscopy or flow cytometry.

2.2 Adherent cells

(1) Culture adherent cells on sterile coverslips.

(2) Remove the coverslip from the culture medium, suck out the excess culture medium, and place the coverslip in a humid environment.

(3) Add 100µL of Hoechst dye working solution from one corner of the coverslip, shake gently to evenly cover all cells with the dye, and incubate at room temperature in the dark for 5-15 minutes.

(4) Aspirate away the dye working solution, wash the coverslip 2 to 3 times with culture solution, and observe with a fluorescence microscope.

Precautions:

- ① For staining of fixed cells or tissue samples, the fixative needs to be removed by rinsing after fixation;
- ② Hoechst 33258 staining is usually performed after other staining. If no other staining is required, Hoechst 33258 staining is performed directly;
- ③ To slow down fluorescence quenching, it is recommended to use anti-fluorescence quenching mounting agent;
- ④ Fluorescent dyes all have quenching problems, so please try to avoid light;
- ⑤ Hoechst 33258 is irritating to the human body. For your safety and health, please wear a lab coat and disposable gloves when operating.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

References:

[1]. K D Harshman, P B Dervan. Molecular recognition of B-DNA by Hoechst 33258. 1985 Jul 11;13(13):4825-35. doi: 10.1093/nar/13.13.4825.

Background

The nucleic acid stain Hoechst 33258 (Ex/Em: 352/461 nm) is frequently utilized as a cell-permeable nuclear counterstain that emits a blue fluorescence upon binding to dsDNA. Hoechst 33258 is commonly employed in various studies related to cell counting, cell cycle analysis, and cell replication. It is particularly useful in identifying condensed nuclei in apoptotic cells, as well as in combination with BrdU staining for cell-cycle studies. Hoechst 33258 works similarly to Hoechst 33342 but is less cell permeable.

Hoechst dyes are also useful for monitoring cell viability by tracking changes in their emission spectra. As minor groove-binding DNA stains with AT selectivity, the Hoechst dyes are able to bind to all nucleic acids, but they show a greater fluorescence enhancement for AT-rich double-stranded DNA strands compared to GC-rich strands [1]. This property has been exploited to identify Q-bands in chromosomes, which are regions rich in AT base pairs that fluoresce brightly when stained with the quinacrine dye [2].

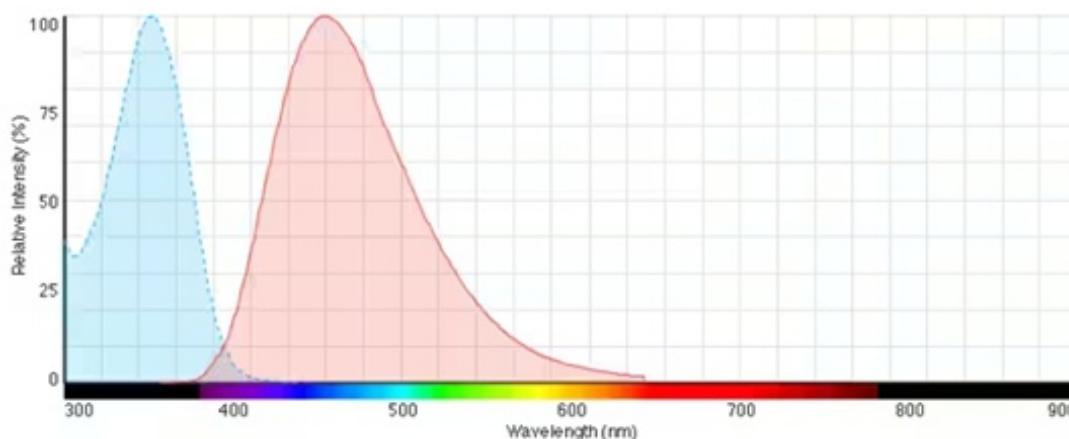


Fig. Fluorescence excitation and emission spectra of Hoechst 33258 bound to DNA

References:

[1]. Portugal J, Waring MJ. Assignment of DNA binding sites for 4', 6-diamidine-2-phenylindole and bisbenzimidazole (Hoechst 33258). A comparative footprinting study. Biochimica et Biophysica Acta (BBA)-Gene Structure and Expression. 1988 Feb 28;949(2):158-68.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

[2]. Weisblum B, Haenssler E. Fluorometric properties of the bibenzimidazole derivative Hoechst 33258, a fluorescent probe specific for AT concentration in chromosomal DNA. Chromosoma. 1974 Sep;46(3):255-60.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA