

Product Data Sheet

Product Name: Ionomycin free acid

Cat. No.: GC15446

Chemical Properties

Cas. No. 56092-81-0

Chemical Name (4R,6R,8R,10E,12S,14S,16E,18S,19S,20R,21S)-11,19,21-trihydroxy-22-((2R,2'R,5R,5'R)-5'-((R)-1-hydroxyethyl)-2,5'-dimethyloctahydro-[2,2'-bifuran]-5-yl)-4,6,8,12,14,18,20-h eptamethyl-9-oxodocosa-10,16-dienoic acid

SMILES O[C@H](C)[C@]1(C)O[C@H](CC1)[C@]2(C)O[C@@H](C[C@@H]([C@@H](C)[C@H]([C@H]/C=C/C[C@H](C)C[C@H](C)/C(O)=C\C([C@H](C)C[C@H](C)C[C@H](C)CCC(O)=O)=O)C)O)O)CC2


Formula C₄₁H₇₂O₉ M.Wt 709.01

Solubility 1.4mg/mL in DMSO, 2.5mg/mL in DMF Storage Desiccate at -20°C

General tips For obtaining a higher solubility , please warm the tube at 37 °C and shake it in the ultrasonic bath for a while. Stock solution can be stored below -20°C for several months.

Shipping Condition Evaluation sample solution : ship with blue ice All other available size: ship with RT , or blue ice upon request.

Structure

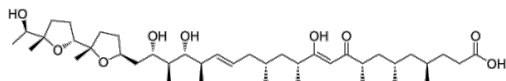

Protocol

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

Cell experiment [1]:	
Cell lines	L6 cells (rat skeletal muscle myoblast cell line)
Preparation Method	L6 cells were maintained in DMEM supplemented with 10% FBS and 5% penicillin/streptomycin at 37°C, 5% CO ₂ . Detached L6 cells were treated with 1µM Ionomycin free acid in PBS with 1.5mM Ca ²⁺ for 10min at room temperature, then rinsed with Ca ²⁺ -free PBS and resuspended in 8.5% sucrose and 0.3% dextrose solution.
Reaction Conditions	1µM; 10min
Applications	Ionomycin free acid significantly increased cytoplasmic conductivity and decreased cytoplasmic permittivity
Animal experiment [2]:	
Animal models	Female athymic nude mice (BALB/c nu/nu, 6-8 weeks old)
Preparation Method	Mice were subcutaneously inoculated with 5×10 ⁶ HT1376 human bladder cancer cells in the flank. Starting 5 days after implantation, mice received intratumoral injection of 100µg Ionomycin free acid in 100µL PBS 3 times per week for 4 weeks. Tumor size (length×width) was measured with caliper at least twice weekly. In a separate combination arm, mice additionally received intraperitoneal injection of 50µg cisplatin in 100µL PBS 3 times per week for 4 weeks. Tumors were resected at 2 and 4 weeks for Northern blot analysis of Bcl-2 and Bax mRNA.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

Dosage form	100µg per mouse; intratumoral; 3 times/week for 4 weeks
Applications	Ionomycin free acid reduced subcutaneous HT1376 tumor size in nude mice. Ionomycin free acid downregulated Bcl-2 mRNA and upregulated Bax mRNA in HT1376 tumor tissues. Combined intratumoral Ionomycin free acid plus intraperitoneal cisplatin further suppressed HT1376 tumor growth compared with either agent alone.
References: [1] Ferguson CA, Santangelo C, Marramiero L, et al. Broadband electrical spectroscopy to distinguish single-cell Ca^{2+} changes due to ionomycin treatment in a skeletal muscle cell line. <i>Sensors</i>. 2023 Apr 28;23(9):4358. [2] Miyake H, Hara I, Yamanaka K, et al. Calcium ionophore, ionomycin inhibits growth of human bladder cancer cells both in vitro and in vivo with alteration of Bcl-2 and Bax expression levels. <i>J Urol</i>. 1999 Sep;162(3):916-21.	

Background

Ionomycin free acid is a highly selective calcium ionophore produced by *Streptomyces globatus*. Ionomycin free acid binds Ca^{2+} and facilitates transport across biological membranes. Ionomycin free acid elevates cytosolic Ca^{2+} levels to activate calcium-dependent signaling pathways, induce apoptosis and autophagy. Ionomycin free acid is often used together with PMA to activate T cells and promote the production of cytokines such as IL-2. Ionomycin free acid can be used for research related to calcium signal transduction, apoptosis, immune activation and cytokine production^[1-4].

In vitro, 1µM Ionomycin free acid treatment for 16 hours reduced degranulation and cytotoxic activity of NK cells against K562, Jurkat, Molt4 and 721.221 target cells. Ionomycin free acid decreased IFN-γ production and impaired the ability of NK cells to form stable conjugates with target cells^[5]. 1µM Ionomycin free acid treatment for 10 minutes in phosphate-buffered saline containing 1.5mM Ca^{2+}

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

elevated cytosolic Ca^{2+} concentration in L6 cells. Ionomycin free acid significantly increased cytoplasmic conductivity and decreased cytoplasmic permittivity^[6]. 5-10 μM Ionomycin free acid treatment for 2-10 minutes rapidly and persistently elevated cytosolic free Ca^{2+} concentration in rat kidney cortical epithelial cells. Ionomycin free acid dissipated mitochondrial membrane potential and caused loss of cell viability^[7].

In vivo, a single injection of 6 μL 25 μM Ionomycin free acid into the peroneal and tibial fascicles of the left sciatic nerve of Sprague-Dawley rats induced vesicular demyelination originating at paranodes and Schmidt-Lanterman incisures and extending along compact internodal myelin^[8]. A single injection of 0.5 μL 2-50 μM Ionomycin free acid at depths of 400 μm and 700 μm at the T10-T11 intervertebral space of Sprague-Dawley rats induced focal primary vesicular demyelinating lesions in the dorsal column, accompanied by oligodendrocyte damage, astrocyte edema and axonal degeneration^[9]. Intratumoral injection of 100 μg Ionomycin free acid 3 times per week for 4 weeks in BALB/c nu/nu mice subcutaneously inoculated with 5×10^6 HT1376 cells reduced subcutaneous HT1376 tumor size. Ionomycin free acid downregulated intratumoral Bcl-2 mRNA expression and upregulated Bax mRNA expression^[10].

References:

- [1] Phillippe M, Chien EM, Freij M, et al. Ionomycin-stimulated phasic myometrial contractions. Am J Physiol. 1995 Oct;269(4 Pt 1):E779-85.
- [2] Chen C, Sun T, Yin M, et al. Ionomycin-induced mouse oocyte activation can disrupt preimplantation embryo development through increased reactive oxygen species reaction and DNA damage. Mol Hum Reprod. 2020 Oct 1;26(10):773-783.
- [3] Aiken EJ, Kilberg BG, Yu S, et al. Ionomycin-Induced Changes in Membrane Potential Alter Electroporation Outcomes in HL-60 Cells. Biophys J. 2018 Jun 19;114(12):2875-2886.
- [4] Hunt JM, Silinsky EM. Ionomycin-induced acetylcholine release and its inhibition by adenosine at frog motor nerve endings. Br J Pharmacol. 1993 Oct;110(2):828-32.
- [5] Romera-Cardenas G, Thomas LM, Lopez-Cobo S, et al. Ionomycin treatment renders NK cells hyporesponsive. PLoS ONE. 2016 Mar 23;11(3):e0150998.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA

Product Data Sheet

- [6] Ferguson CA, Santangelo C, Marramiero L, et al. Broadband electrical spectroscopy to distinguish single-cell Ca^{2+} changes due to ionomycin treatment in a skeletal muscle cell line. *Sensors*. 2023 Apr 28;23(9):4358.
- [7] Jiang T, Grant RL, Acosta D. A digitized fluorescence imaging study of intracellular free calcium, mitochondrial integrity and cytotoxicity in rat renal cells exposed to ionomycin, a calcium ionophore. *Toxicology*. 1993;85(1):41-65.
- [8] Smith KJ, Hall SM. Peripheral demyelination and remyelination initiated by the calcium-selective ionophore ionomycin: in vivo observations. *J Neurol Sci*. 1988 Jun;83(1):37-53.
- [9] Smith KJ, Hall SM. Central demyelination induced in vivo by the calcium ionophore ionomycin. *Brain*. 1994;117(6):1351-6.
- [10] Miyake H, Hara I, Yamanaka K, et al. Calcium ionophore, ionomycin inhibits growth of human bladder cancer cells both in vitro and in vivo with alteration of Bcl-2 and Bax expression levels. *J Urol*. 1999 Sep;162(3):916-21.

Caution: Product has not been fully validated for medical applications. For research use only.

Tel: (909) 407-4943 Fax: (626) 353-8530 E-mail: tech@glpbio.com

Address: 10292 Central Ave. #205, Montclair, CA, USA